


Maxi Sky 2

Kit N° 700-14464 for "G" Track Profile

Parts included:


MS2 Lift for "G" Track Profile	Maxi Sky 2 # MSX11-GX-XX
Charging Station for Other Tracks	# 700-15522-1
Guiding-wheels assembly for "G" Track Profile	# 700.14464
Universal End Stopper	# A5320-C

1. Install the guiding-wheels as shown in the picture below.

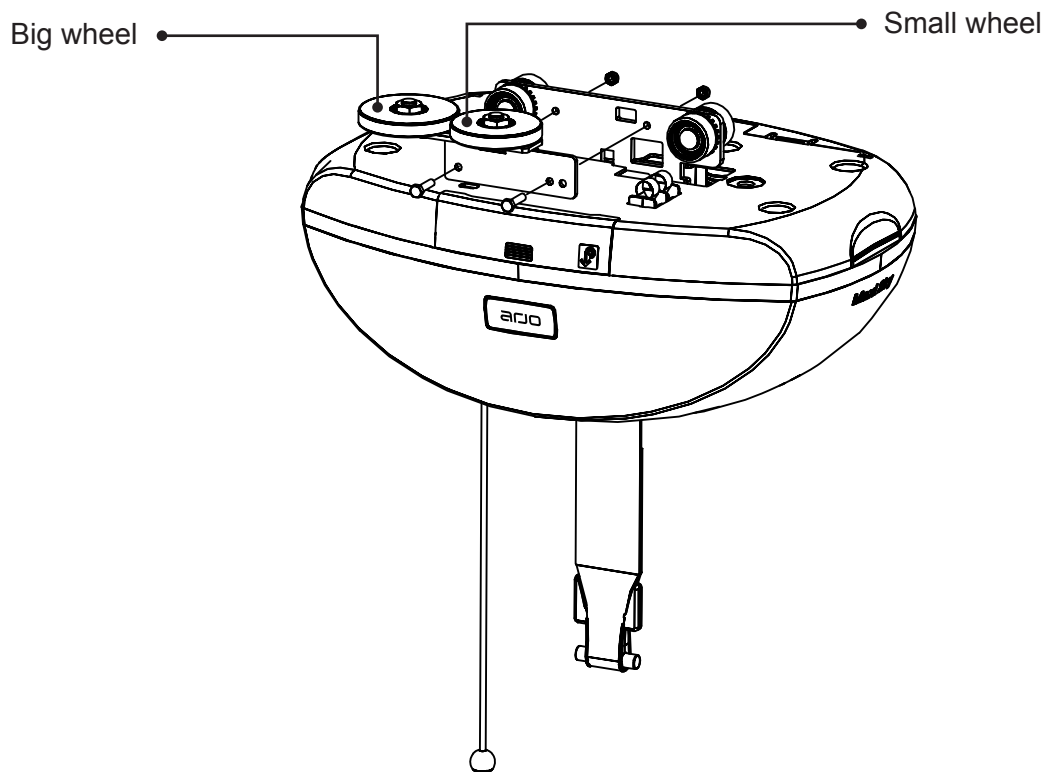

CAUTION

Make sure that the wheels are centered with the trolley and the lift.

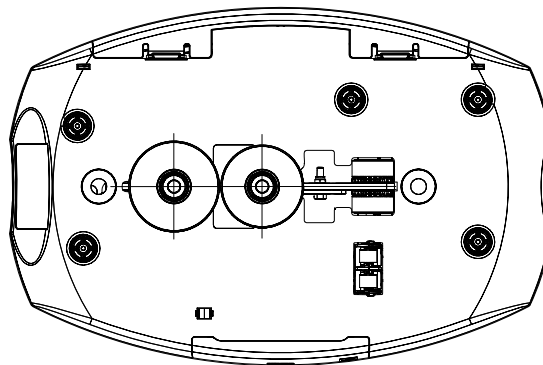

2. Install the end stopper, Charger and Maxi Sky 2 according to the installation procedures provided.