


Sara Plus wipeable standing sling



Wipe to
clean

With no need to launder, and designed to be wiped clean between patients, Arjo's wipeable standing c-hook sling is easy to use, easy to clean, and minimises downtime between uses.

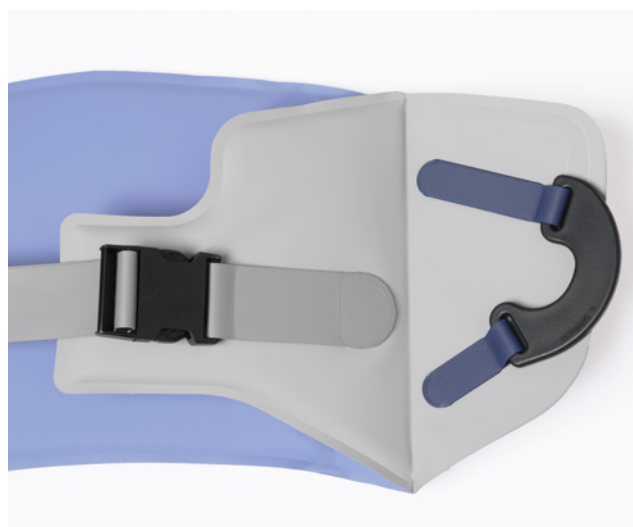

Design benefits

Increased operational efficiency

An excellent alternative for healthcare facilities that prefer to avoid the laundry process and the time required for its management.

With this wipeable sling, the caregivers can transfer the sling from one patient to another with minimum delay, able to clean and disinfect the sling just by wiping it down between uses. The sling is always accessible with the lift.

Prevent cross-contamination

With no stitches, the product is manufactured with welded seams only and the body, the belt and the c-hooks of the sling are made of non-porous material, which makes it efficient to clean and decontaminate with a wipe, greatly reducing the risk of cross contamination.

Increase cost effectiveness

No more lost or damaged slings resulting the laundering process.

Compliance to disinfection protocol also helps to reduce costs related to cross contamination and healthcare-associated infections (HAI).

Key features

SECURE COMFORTABLE FIT

The sling material combined with the motion of Sara Plus during the sit to stand process is designed to keep the sling in place during transfer to support patient comfort.

COLOUR CODED SIZING

All slings are size colour-coded for ease when choosing the right sling.

NON POROUS MATERIAL

The body, the belt and the c-hooks of the sling are made of non-porous material and hence completely wipeable.

WELDED SEAMS

Non-stitched, welded seams prevent potential fluid ingress and contamination of the sling.

C HOOK

For simple integration with the Sara Plus

FIXED BELT


The belt and sling are permanently fixed to prevent the two being separated and the belt being lost.

CONCAVE SHAPE

The concave shape at the back of the sling enables access to the patient for undressing and hygiene care.

SMOOTH PATIENT SURFACES

The inside of the sling is free from seams, straps and other components, creating a smooth surface towards the patient's body.

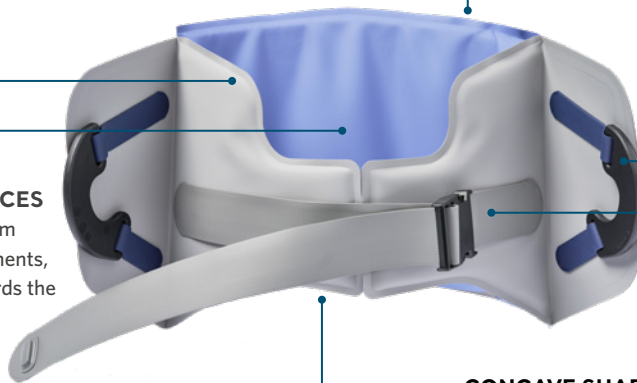

PRODUCT SPECIFICATIONS

General

Compatible lifts	Sara® Plus
Safe Working Load, SWL	200 kg / 440 lb

Size (Article Number)	Waist size in cm (in)
Small (KKA6000-S)	64-84 cm (25" - 33")
Medium (KKA6000-M)	76-104 cm (30"- 41")
Large (KKA6000-L)	88-124 cm (34 ½" - 49")
Extra Large (KKA6000-XL)	108-160 cm (42 ½"- 63")

Cleaning/disinfection

This sling can be wiped down with the following types of cleaning/disinfection agents:	Quaternary ammonium, Isopropyl alcohol, Phenolic, Bleach (Sodium hypochlorite)*, Peracetic acid, Hydrogen peroxide.
Please refer to the "Instructions for use" for more information.	

**Sling may show some signs of discolouration and become slightly sticky after repeated contact with bleach.*