

HEIGHT-ADJUSTABLE RECLINING

Parker



Parker® — a multipurpose hygiene experience. The wide side-entry door makes ingress and egress easy for individuals with limited lower-body joint mobility. *Parker* can be used as a relaxing air spa or immersion bath, a shower cabinet and/or a podiatry tank for foot care — all with person-centered care in mind.

Designed with Albert, Barbara and Carl in mind. Also suitable for Doris when paired with a passive lift.



Smooth and ergonomic transfers

Parker features a side door that swings up to facilitate unprecedented access for more active residents. *Parker* can also be used with a range of active or passive ceiling lifts, depending on the mobility level of the individual resident. This enables smooth transfers without the need for manual lifting for a wide range of resident mobility levels.

Anti-scalding protection features

With people in mind, *Parker* includes clear temperature displays, temperature- and pressure-regulated mixing valves, and automatic hot water shut-off when exceeding safe temperatures.

Wellness features

Select from a still immersion bath, Air Spa™ and/or Sound & Vision™ to enable you to create a calm, multisensory environment to help your residents and residents relax.

Key features

WIDE SIDE

The wide door opening is designed to integrate with a range of transfer aids, allowing for easy, comfortable access for active or semi-active residents.

HEIGHT-ADJUSTABLE AND RECLINING

Allows caregiver to reach preferred working height without compromising resident comfort.

ANTI-SCALDING PROTECTION FEATURES

Water temperature displays, temperature- and pressure-regulated mixing valves, and solenoid shut-off for unsafe temperatures provide the safest environment for the resident.

SHOWER WAND WITH LOCKABLE TRIGGER

Enables residents to participate with their bath and showering activities.

HEAD AND BACK CUSHIONS

For additional comfort.

INTEGRATED CLEANING AND DISINFECTION SYSTEM

AIR SPA (OPTIONAL)

With nine speed-controlled air jets.

AUTO-FILL

Auto-fill provides efficiency — allowing the tub to be filled while the resident is retrieved for the bath. It aids in providing a quiet, calm bathing experience.

SOUND & VISION (OPTIONAL)

For an enhanced multisensory experience.



GENERAL

Maximum safe working load (SWL) (patient + water)	661 lbs (300 kg)
Maximum patient weight	463 lbs (210 kg)
Total length in reclined position	77 1/4" (1,960 mm)
Total width	29 1/2" (750 mm)
Height of door with bathtub in highest position, feet in lowest position	98 7/8" (2,510 mm)
Seat height at lowest position	22 1/4" (565 mm)
Water consumption (prefilled foot well volume/auto-fill level)	25.4 gallons (96 liters)
Water consumption up to overflow level	60.8 gallons (230 liters)
Filling time at 3 bar (43.5 PSI) dynamic pressure (prefill)	1.7 minutes
Emptying time from auto-fill level	3 minutes

STANDARD FEATURES


- Bath and shower temperature display
- Temperature- and pressure-regulated mixing valve
- Automatic hot water shut-off when exceeding safe temperature
- Shower wand with lockable trigger
- Auto-fill
- Overflow
- Tamper-proof security lock for control panel functions
- Battery backup for emergency tilting and lowering
- Emergency stop
- Head support cushion

OPTIONAL FEATURES

- Footrest and back/seat support cushions
- Air Spa therapy
- Sound & Vision therapy
- Tub accent color available in blue or gray

SPACE REQUIREMENT

Blue area denotes the required working area to facilitate transfer and hygiene activities for residents with mobility levels A to C while providing adequate access for resident, equipment and assisting caregiver. Additional space might be required to facilitate transfer of resident with mobility level D using passive lift.

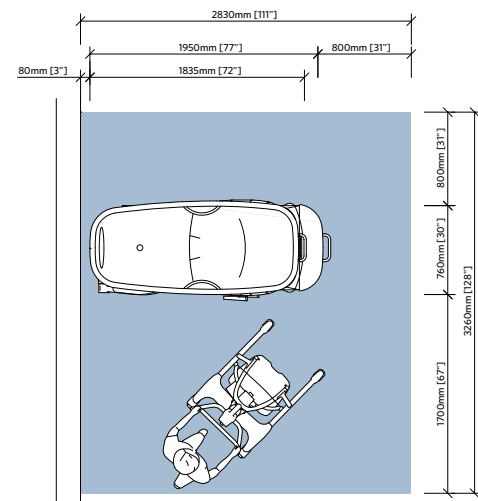

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure ulcers and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind. www.arjo.com

Arjo Inc. · 2349 West Lake Street · Addison, Illinois 60101 · 800-323-1245 · www.arjo.us

Arjo Canada Inc. · 90 Matheson Blvd West, Suite 350 · Mississauga, ON, L5R 3R3 · 800-665-4831 · Info.Canada@arjo.com · www.arjo.ca

arjo