

# Minstrel



Minstrel is the classic loop sling lifter for everyday transfer routines with dependent residents.

## ARJO MOBILITY GALLERY



D



E

© Arjo 2019

D. Unable to bear weight through feet. Is able to sit if well supported.  
E. Difficulty maintaining sitting, hence might remain lying/partially lying in bed, and may only sit in a special supportive chair.  
Please contact Arjo for further information.

- **Simplicity**  
Easy to understand and simple to operate.
- **Eye contact**  
The handset enables the caregiver to move and stay close to the resident and maintain good contact throughout the lift and transfer process.
- **Loop slings**  
Loop sling solutions are available in a wide range of designs and sizes.
- **Lifting from floor**  
With Minstrel, a single caregiver can move a resident who has fallen on the floor.
- **Electronic scale (option)**  
Minstrel is available with an optional scale that enables the caregiver to weigh the resident while performing a transfer.

# Key features

## SOLID HANGER BAR

Solid hanger bar assembly designed to minimise swinging of the resident during the lift and transfer.

## LOW-FRICTION CASTORS

Castors designed to roll easily and enable the caregiver to manoeuvre during the transfer.

## SCALE (OPTION)

Optional suspended scale supports an efficient weighing routine for dependent residents.

## DUAL CONTROLS

Minstrel can be operated from the control panel on the mast or with the hand control, enabling the caregiver to stay close to the resident while operating the lifter.

## CHASSIS ADJUSTMENT

The chassis can be adjusted to allow access and manoeuvrability around large chairs, wheelchairs, beds and toilets.



## PRODUCT INFORMATION

Max safe working load	190 kg	419 lb
Lifter weight (with spreader bar and battery)	56 kg	123 lb
Max lifting height	2000 mm	78 ¾ in
Min lifting height	625 mm	24 ½ in
Lifting stroke	1375 mm	54 ⅞ in
Min total storage height	1330 mm	52 ⅞ in
External width legs closed	645 mm	25 ⅜ in
External width legs open	1181 mm	46 ½ in
Internal width leg closed	545 mm	21 ½ in
Internal width leg open	1080 mm	42 ½ in
Leg clearance (floor to top of legs)	110 mm	4 ⅓ in
Chassis clearance (floor to bottom of chassis)	45 mm	1 ¾ in
Lifter length	1275 mm	50 ¼ in
Turning diameter	1420 mm	56 in
Lifter protection class	IPX4	
Hand control protection class	IPX7	
Battery	24V, 4 Ah	
Battery charge indicator		
Service meter – displays the total time in use (hours)		
Emergency stop and system failure override		
Automatic safety cut-out if lowered onto an obstacle		
Low-friction castors, the rear two with brakes		


## ACCESSORIES

Wide range of special purpose loop slings and flites (disposable slings)

Option: suspended scale

Arjo disinfectants

## SPACE REQUIREMENT

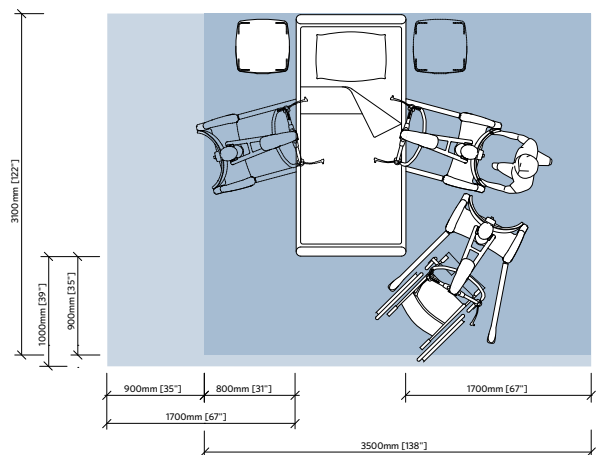

- Blue area shows the minimum working area required to enable the caregiver to use the floor lifter ergonomically from one side.
- Light blue area shows extension of working area required to facilitate activities from either side to provide adequate access for the resident, floor lifter and caregiver.