Miranti[®]



Fully compatible with a range of Arjo recumbent baths, Miranti is designed for the needs of patients and residents with the mobility level of Doris or Emma, who need to remain in a lying position. Miranti can be operated by a single caregiver while allowing them to maintain ergonomically sound working positions throughout the task.

Designed with Doris and Emma in mind





Enables recumbent residents to enjoy the benefits of bathing

Featuring a full-length stretcher that provides full body support, Miranti can be used to collect a resident from their bedside or wheelchair and transfer them directly into the tub, allowing residents who need to remain in a lying position to enjoy the benefits of bathing.

Designed to be operated by a single caregiver

Miranti is height-adjustable and features a fifth wheel power drive, central braking, and a powered back rest which can be raised or lowered to the resident's preference.

Designed for easy cleaning and disinfection

All external surfaces of the Miranti are designed to be wipeable, including the stretcher which is coated in moulded, high-density polyurethane foam.



Key features

HEIGHT-ADJUSTABLE

With wide stroke range to reach preferred height for both caregiver and resident.

ERGO GRIPS

Situated at both ends of the stretcher for manoeuvring and power drive control.

FOLDABLE GRIP & SIDEGUARD

Provides additional security and support, and folds away to facilitate transfers.

POWERED BACK REST

Can be raised or lowered at the touch of a button.

DUAL CONTROLS

Enables caregiver to operate the device while staying close to resident.



ADJUSTABLE NECK AND THIGH SUPPORT CUSHIONS

REVERSIBLE STRETCHER

Either end of the stretcher can be used as the head-end to enable access to either side of the bed or bath.

SCALE (OPTIONAL)


Class III integrated electronic scale offers a convenient alternative for weighing residents.

FIFTH WHEEL POWER DRIVE (OPTIONAL)

Helps to reduce workload when transferring and manoeuvring.

CENTRAL BRAKE

SPACE REQUIREMENT

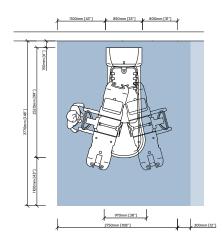

Required working area to facilitate bathing and transfer activities from either side of a long recumbent tub while providing adequate access for the resident, equipment and assisting carer.

STANDARD FEATURES & ACCESSORIES

Height-adjustable

Powered back rest

Central brakes

Dual controls on chassis and handset

Emergency raising and lowering

Foldable security grip and sideguard

Ergo grips at head and foot ends

Reversible stretcher

Wipeable safety belt

Adjustable neck and thigh support cushions

Battery charger and 2x rechargeable batteries

PROD	UCT	INFOR	MATION	

Safe working load	160 kg (353 lb)
Total weight (lift, power drive and scale)	103 kg (227 lb)
Total width	930 mm (36 % in)
Total length	1970 mm (77 ½ in)
Stretcher height, lowest position	475 mm (18 ¾ in)
Stretcher height, highest position	1065 mm (41 % in)
Lifting stroke range (stretcher)	590 mm
Height of lift, max (back rest up)	1890 mm (74 ¾ in)
Height of lift, min (back rest up)	1300 mm (51 ½ in)
Chassis width, max	950 mm (37 % in)
Chassis width between legs, min	382 mm (15 in)
Clearance under chassis	30 mm (1 ¼ in)
Chassis length, max (incl. brakes)	905 mm (35 % in)
Turning radius (flat stretcher position)	ø 1990 mm (78 ½ in)
Supply voltage	24 V DC
Max power	150 VA
Protection class	IP X4
Protection class (handset)	IP X7

OPTIONAL FEATURES

Integrated weighing scale

Fifth wheel power drive

COMPATIBLE BATH TUBS

Rhapsody standard (System 23) and long (System 25) tub lengths

Primo standard (System 23) tub length

Classic Line

For more detailed information visit www.arjo.com

