



Designed for
the miracle
of life

MATERNITY CARE RANGE

Prioma[®] Birthing Bed



Reliability and comfort for one of life's greatest moments



The Prioma Birthing Bed is an electrically operated delivery bed designed to provide comfort and support to the mother, as well as efficient and reliable functionality for the delivery team throughout all stages of the birthing process.

The soft appeal in the design, functionality and ease-of-use provide a welcoming, practical and affordable solution for any labour suite or maternity ward.



Converts easily from a standard bed to a versatile delivery platform.



Standard features

- Safe Working Load (SWL) 200 Kg
- 150 mm twin central locking castors
- Linak actuators
- Easily removable ABS backrest panel
- Easily removable high pressure compact laminate seat & foot panels
- Head end drop down side rails
- Embedded controls for nurse and patient
- Integral bed extension - 580 mm
- Battery back up
- Retracting backrest
- Manual CPR backrest quick release
- Removable lightweight head & foot panels
- Trendelenburg tilt 14°, reverse Trendelenburg tilt 12°

Optional features

- Under bed light - set in permanent on-mode
- 125 mm twin central locking castors
- Linak nurse control handset

A versatile, comfortable and functional choice



Fully integrated patient and nurse control panels



Large under tray



Ergonomically designed, sturdy and reliable leg stirrups



Retracting backrest

Standard accessories

A range of attachable accessories come as standard to optimise positioning for delivery and clinical procedures.

- Under tray
- Grab handles
- Footrest plate
- Upper and lower leg supports
- Leg stirrups
- Mattress sections



Lower foot supports



Grab handles



Footrest plate

Optional accessories

- Storage trolley - powder coated finish
- Self help pole - includes strap and handle, available in stainless steel or powder coated finish
- IV pole - stainless steel finish
- "C" size oxygen bottle holder - clamps to head end of bed
- Extra set of mattress covers
- Dust cover




Optional storage trolley

Optimal positioning to aid in comfort and support clinical procedures



Bed positions

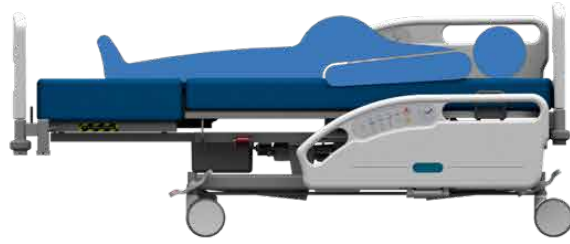



Basic Dimensions

Overall length - bed configuration	2130 mm
Overall length - obstetric configuration	1580 mm
Bed extension length	580 mm
Mattress surface length - between panel faces	2010 mm
Overall width	1000 mm
Deck height range - 150mm twin central locking castors	485 mm to 840 mm +/- 5 mm
Trendelenburg and reverse Trendelenburg tilt - electronic	0 to 14° (12° reverse Trendelenburg)
Backrest angle adjustment - electronic	0 to 70°
Mattress size seat and backrest section	850 mm (w) x 1350 mm (l) x 125 mm (h)
Mattress size extension bolster	850 mm (w) x 750 mm (l) x 170 mm (h)

Product Weight (approximate)

Bed without accessories	160 Kg
Mattress sections	14 Kg
Safe Working Load (SWL)	200 Kg
Maximum patient weight	135 Kg

Electrical Data

Power in - 100 - 240V ~ 50/60Hz Max 3.9A
Duty rating - intermittent 10%, max 2 min / 18 min
Electric Shock Protection - Class 1 Type B
Liquid ingress protection - IPX4
Internally powered equipment - battery backup only - 2 x 12V series connected batteries, sealed chargeable lead/Acid Gel 1.2 amp hours

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure ulcers and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind. www.arjo.com

Arjo AB, Hans Michelsensgatan 10, SE-211 20 Malmö, Sweden. Phone: +46 (0) 10 335 45 00

www.arjo.com