

ACCESSORY BROCHURE

Minuet 2

Community bed range





Accessories to help you provide better care

The Minuet 2 community bed range has an extensive offering of accessory items, which both complement the bed and provide additional benefits for residents and carers. Designed specifically for the Minuet 2 community bed range, these accessories reflect the same high quality and standards as the beds themselves.

Minuet 2 accessories can be ordered independently of the bed, and can be fitted retroactively to the bed as required (unless otherwise stated). All accessories have been fully tested to ensure compatibility with the Minuet 2 bed range.

The Minuet 2 community bed range has been independently tested by the globally renowned organisation UL, to certify that it is compliant with the latest stability, dimensional and safety requirements of the 2-52 Standard.



Summary

| ACCESSORY TYPE | ACCESSORY ITEM | PART NUMBER | PAGE |
|----------------------------------|---|---|----------------------------|
| Repositioning and egress support | Egress assist rail, foldable, standard height Egress assist rail, foldable, extra height Mobility support handle (left) Mobility support handle (right) | CM-ACC22-1.52 CM-ACC23-1.52 CM-ACC11-1 CM-ACC12-1 | 4 4 4 4 |
| Wooden safety sides | Wooden, full length safety sides, standard height, (grey/beech)* Wooden, full length safety sides, extra height, (brown/beech)* Wooden, full length safety sides, standard height (grey/beech) Wooden, full length safety sides, standard height (brown/beech) Wooden, full length safety sides, standard height (grey/cherry) Wooden, full length safety sides, extra height (grey/beech) Wooden, full length safety sides, extra height (brown/beech) | CM-ACC00-1 CM-ACC00-3 CM-ACC00-1.52 CM-ACC00-3.52 CM-ACC01-1.52 CM-ACC00-1.14 CM-ACC00-3.14 | 5 5 5 5 5 6 |
| Steel safety sides | Steel, full length safety sides, standard height Steel, 3/4 length safety sides, standard height, folding Steel, 3/4 length safety sides, extra height folding | CM-ACC04-1.52 CM-ACC06-1.52 CM-ACC07-1.52 | 7 7 7 |
| Safety sides height extensions | Height extensions for wooden full length safety sides (grey) Height extensions for wooden full length safety sides (brown) | CM-ACC02-1.52 CM-ACC02-3.52 | 8 |
| Safety sides pads | Pads for full length safety sides Pads with top cushion for full length safety pads | CM-ACC03-1 CM-ACC26-1 | 9 |
| Lifting poles | Lifting Pole with strap and handle, single position | CM-ACC09-1 | 10 |
| IV pole | IV pole | CM-ACC14-1 | 10 |
| Other | Mattress Pump Bracket Urine Bag Holder Mattress Retainers (supplied with the bed) Flexible handset holder | CM-ACC15-1 CM-ACC17-1 CM-ACC24 CM-ACC13-1 | 11 11 11 11 |

^{*} Provided for backwards compatibility with beds sold before 2013 that are compliant with the older IEC 60601-2-38 standard for safety of medical beds. These cannot be sold with new beds compliant with the succeeding IEC 60601-2-52 standard.

Repositioning and egress supports



CM-ACC22-1.52

Egress assist handle, foldable, standard height

In the raised position, it can be used to aid a resident in rising from the bed. When lowered, it provides additional security for residents in cases where the use of longer side rails is not desirable. Supplied as a pair.

- Product weight is 6 kg (13.3 lb) each
- Product length is 745 mm (29.3 in)
- Product height above MSP is 400 mm (15.75 in)
- Compatible with all Minuet 2 beds



CM-ACC23-1.52

Egress assist handle, foldable, extra height

Function is the same as for the standard height foldable egress rail CM-ACC22-1.52.

Extra height egress assist rails should be used when the mattress thickness is between 150 mm and 172 mm. Supplied as a pair.

- Product weight is 7.3 kg (16 lb) each
- Product length is 785 mm (30.9 in)
- Product height above MSP is 510 mm (20 in)
- Compatible with all Minuet 2 beds



CM-ACC11 (left) / CM-ACC12 (right)

Mobility support handle

The mobility support handle is used to assist the resident in getting on or off the bed. It is fitted in the accessory socket at the head end of the bed.

CM-ACC11 fits on the left side of the bed (when viewed from the head end) and CM-ACC12 fits on the right side of the bed.

- Safe working load is 75 kg (165 lb)
- Product weight is 4 kg (8.8 lb)
- Compatible with all Minuet 2 beds

Wooden safety sides



CM-ACC00-1

Wooden full length, standard height, safety sides (grey/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Standard height safety sides should be used when the mattress thickness does not exceed 150 mm.

- Product weight is 10 kg per side (22 lb)
- Product height above MSP is 410 mm (16 in)
- Compatible with 2-38 compliant beds only*



CM-ACC00-3

Wooden full length, extra height, safety sides (brown/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Extra height safety sides should be used when the mattress thickness is between 150 mm and 172 mm.

- Product weight is 10 kg per side (22 lb)
- Product height above MSP is 410 mm (16 in)
- Compatible with 2-38 compliant beds only*



CM-ACC00-1.52

Wooden full length, standard height, safety sides (grey/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Standard height safety sides should be used when the mattress thickness does not exceed 150mm.

- Product weight is 10 kg per side (22 lb)
- Product height above MSP is 410 mm (16 in)
- Compatible with all Minuet 2 beds



CM-ACC00-3.52

Wooden, full length safety sides, standard height, (brown/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Standard height safety sides should be used when the mattress thickness does not exceed 150 mm.

- Product weight is 10 kg per side (22 lb)
- Product height above 410 mm (16 in)
- Compatible with all Minuet 2 beds



CM-ACC01-1.52

Wooden, full length safety sides, standard height, (grey/cherry)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Standard height safety sides should be used when the mattress thickness does not exceed 150 mm.

- Product weight is 10 kg per side (22 lb)
- Product height above MSP is 375 mm (14.8 in)
- Compatible with all Minuet 2 beds



CM-ACC00-1.14

Wooden, full length safety sides, extra height, (grey/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Extra height safety sides should be used when the mattress thickness is between 150mm and 172 mm.

- Product weight is 11 kg per side (24 lb)
- Product height above MSP is 410 mm (16 in)
- Compatible with all Minuet 2 beds

| | 2-52 compliant beds only (2013→) ^{1,2} | | 2-38 compliant beds only (→2013) ^{2,3} | | |
|-----------------|---|--------------------|---|--------------------|--|
| | mattress | thickness | mattress thickness | | |
| ACCESSORY ITEM: | 125-150 mm (5-6'') | 150-172 mm (6-7'') | 125-150 mm (5-6'') | 150-172 mm (6-7'') | |
| CM-ACC00-1 | | | • | | |
| CM-ACC00-3 | | | | • | |
| CM-ACC00-1.52 | • | | • | | |
| CM-ACC00-3.52 | • | | • | | |
| CM-ACC01-1.52 | • | | • | | |
| CM-ACC00-1.14 | | • | | • | |

- All Minuet 2 beds are also compliant with the FDA 7-zone entrapment guidelines if used with the proper mattress as indicated in the table above.
- For more information about how to identify a bed as 2-38 or 2-52 compliant, please see page 12 in this document
- 2-38 is an outdated standard and siderails for 2-38 compliant beds are only provided for beds sold before 2013

Wooden safety sides (continued)



CM-ACC00-3.14

Wooden, full length safety sides, extra height (brown/beech)

Full length sides are supplied as identical pairs. They can be attached to either side of the bed.

Extra height safety sides should be used when the mattress thickness is between 125 mm and 172 mm.

- Product weight is 11 kg per side (24 lb)
- Product height above MSP is 410 mm (16 in)
- Compatible with all Minuet 2 beds

Steel safety sides



CM-ACC04-1.52

Steel, full length safety sides, standard height

The full length steel safety sides are supplied as a pair to fit the right and left hand sides of the bed, as viewed from the head end.

The standard height steel full length safety sides can be used when the mattress thickness is between 125 mm and 150 mm.

- Product weight is 12 kg (26.7 lb) per side
- Product height above MSP is 380 mm (14.97 in)
- Compatible with all Minuet 2 beds



CM-ACC06-1.52

Steel, 3/4 length safety sides, folding, standard height

The folding safety sides are supplied as a pair to fit the right and left hand sides of the bed, as viewed from the head end.

The standard height safety sides can be used when the mattress thickness is between 125 mm and

- Product weight is 7.7 kg (17 lb) per side
- Product height above MSP is 380 mm (14.97 in)
- Compatible with all Minuet 2 beds



CM-ACC07-1.52

Steel, 3/4 length safety sides, folding, extra height

The folding safety sides are supplied as a pair to fit the right and left hand sides of the bed, as viewed from the head end.

Extra height safety sides can be used when the mattress thickness is between 150 mm and 172 mm.

- Product weight is 9 kg (19.8 lb) per side
- Product height above MSP is 490 mm (19.3 in)
- Compatible with all Minuet 2 beds

| | 2-52 compliant be | 2-52 compliant beds only (2013→) ^{1,2} | | 2-38 compliant beds only (→2013) ^{1,2,3} | | |
|-----------------|--------------------|---|--------------------|---|--|--|
| | mattress | mattress thickness | | thickness | | |
| ACCESSORY ITEM: | 125-150 mm (5-6'') | 150-172 mm (6-7") | 125-150 mm (5-6'') | 150-172 mm (6-7'') | | |
| CM-ACC04-1.52 | • | | • | | | |
| CM-ACC06-1.52 | • | | • | | | |
| CM-ACC07-1 52 | | • | | • | | |

- All Minuet 2 beds are also compliant with the FDA 7-zone entrapment guidelines if used with the proper mattress as indicated in the table above.
- For more information about how to identify a bed as 2-38 or 2-52 compliant, please see page 12 in this document
- 2-38 is an outdated standard and siderails for 2-38 compliant beds are only provided for beds sold before 2013

Height extensions for safety sides



CM-ACC02-1.52

Height extensions for wooden full length safety sides (grey)

Clamp-on extensions that can be added to the full-length wooden safety sides to increase their height by 140 mm and allow use of a mattress or a combined mattress and overlay combination of a thickness up to 265 mm.

The safety side height extensions are supplied as identical pairs. They can be used on either side of the bed.

- Product weight is 3 kg (6.6 lb) each
- Product height above top of side rail on which it is mounted is 140 mm (5.5 in)
- Compatible with all Minuet 2 full length standard height wooden safety sides



CM-ACC02-3.52

Height extensions for wooden full length safety sides (brown)

Clamp-on extensions that can be added to the full-length wooden safety sides to increase their height by 140 mm and allow use of a mattress or a mattress and overlay combination of a thickness up to 265 mm.

The safety side height extensions are supplied as identical pairs. They can be used on either side of the bed.

- Product weight is 3 kg (6.6 lb) each
- Product height above top of side rail on which it is mounted is 140 mm (5.5 in)
- Compatible with all Minuet 2 full length standard height wooden safety sides

| ACCESSORY ITEM: | CM ACC00-1 | CM-ACC00-3 | CM-ACC00-1.52 | CM-ACC00-3.52 | CM-ACC01-1.52 | CM-ACC00-1.14 | CM-ACC00-3.14 |
|-----------------|------------|------------|---------------|---------------|---------------|---------------|---------------|
| CM-ACC02-1.52 | • | • | • | • | • | | |
| CM-ACC2-3 52 | • | • | • | • | • | | |



CM-ACC03-1

Pads for full length safety sides

These washable pads may be used with any of the full-length safety sides.

They may also be used when height extensions are fitted to the safety sides.

The pads are easy to fit, and can be removed from the bed quickly if the safety sides need to be lowered in an emergency.

- Product weight is 3.9 kg (8.6 lb) per pair
- Compatible with all Minuet 2 full length side rails.



CM-ACC26-1

Pads with top cushion for full length safety sides

These washable pads may be used with any of the full-length safety sides.

They may also be used when height extensions are fitted to the safety sides.

The pads are easy to fit, and can be removed from the bed quickly if the safety sides need to be lowered in an emergency.

Model CM-ACC26-1 provides additional protective padding along the top rail.

- Product weight is 3.9 kg (8.6 lb) per pair
- Compatible with all Minuet 2 standard height full length side rails as well as the full length steel safety sides.

| ACCESSORY ITEM: | CM-ACC00-1 | CM-ACC00-3 | CM-ACC00-1.52 | CM-ACC00-3.52 | CM-ACC01-1.52 | CM-ACC00-1.14 | CM-ACC00-3.14 |
|-----------------|------------|------------|---------------|---------------|---------------|---------------|---------------|
| CM-ACC03-1 | • | • | • | • | • | • | • |
| CM-ACC26-1 | • | • | • | • | • | | |

| ACCESSORY ITEM: | CM-ACC04-1.52 | CM-ACC06-1.52 | CM-ACC07-1.52 | CM-ACC02-1.52 | CM-ACC02-3.52 |
|-----------------|---------------|---------------|---------------|---------------|---------------|
| CM-ACC03-1 | • | | | | • |
| CM-ACC26-1 | • | | | | • |

Lifting pole and IV pole



CM-ACC09-1

Lifting Pole with strap and handle, single position

This fixed lifting pole adopts a centralised position over the bed and assists residents in moving or turning whilst lying on the bed themselves. Height adjustable strap and handle are included.

- Safe working load is 75 kg (165 lb)
- Product weight is 6.7 kg (14.8 lb)
- Compatible with all Minuet 2 beds



CM-ACC14-1

IV note

A straight fixed height IV pole with a single fixed hook. It can be used for drips, pumps and other IV equipment.

- Safe working load is 5 kg (11 lb)
- Product weight is 2 kg (4.7 lb)
- Attaches to accessory socket
- Compatible with all Minuet 2 beds

Other



CM-ACC13-1

Flexible handset holder

The flexible handset holder is designed to position the control handset within easy reach of the resident.

It also allows the controls to be operated by a resident with limited gripping ability.

- Product weight is 1.1 kg (2.7 lbs)
- Compatible with all Minuet 2 beds



CM-ACC17-1

Urine bag holder

The urine bag holder is designed to support a single urine or drainage bag.

It can be fitted to either the left or right-hand side of the bed as required.

- Product weight is 0.1 kg (0.2 lb)
- Compatible with all Minuet 2 beds



CM-ACC15-1

Mattress pump bracket

The pump bracket is designed to support the pump from any compatible Arjo air mattress system.

- Product weight is 0.5 kg (1.1lb)
- Compatible with all Minuet 2 beds




CM-ACC24

Mattress retainers

This accessory comprises four plastic mouldings that are attached to the backrest and leg sections to prevent sideways movement of the mattress.

- Product weight is 0.2 kg (0.44 lb), set of four
- Compatible with all Minuet 2 beds

How to identify 2-52 compliant bed

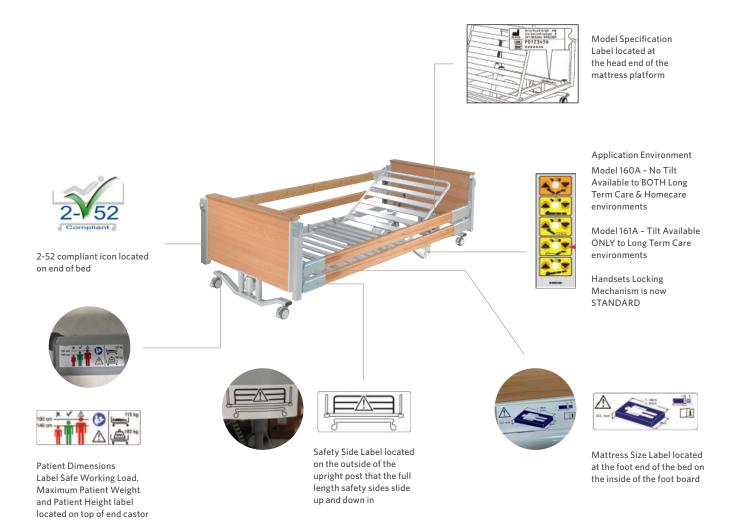

At Arjo, we believe that empowering movement within healthcare environments is essential to quality care. Our products and solutions are designed to promote a safe and dignified experience through patient handling, medical beds, personal hygiene, disinfection, diagnostics, and the prevention of pressure injuries and venous thromboembolism. With over 6500 people worldwide and 65 years caring for patients and healthcare professionals, we are committed to driving healthier outcomes for people facing mobility challenges.

Arjo AB • Hans Michelsensgatan 10 • 211 20 Malmö • Sweden • +46 10 335 4500

www.arjo.com

